

# iSense®

## ELECTRICAL CARDIOMETRY™ SENSORS




# Window to the Circulation®

iSense® sensors are intended for use with the Electrical Cardiometry™ (EC™) monitor AESCULON® and cardiometers ICON® and ICONCore®. iSense® Sensors are 100% disposable and are available in various sizes, for neonate, pediatric, and adults. All sensors utilize a highly conductive AgCl wet gel and strong hypoallergenic adhesion, which provides instant electrical contact for high quality signal conductivity. Four iSense® sensors are packed together in one pouch for the convenience of the user and more importantly keeping the gel fresh for best conductivity and signal quality, especially under difficult patient conditions.



iSense® sensor box

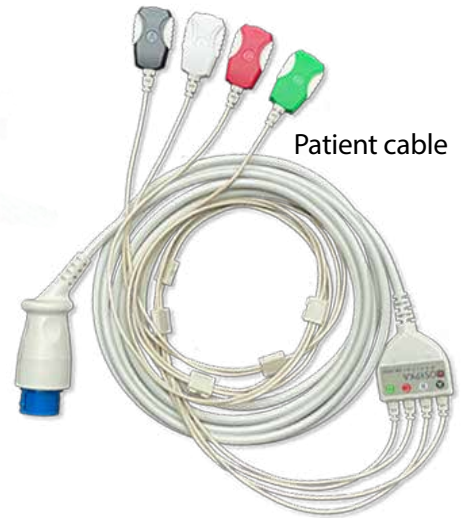

Patient cable

## iSense® Specification:

Sensor	Patient Group	Type/ Size	Article No.	Packaging Unit
	Premature and Neonates	Prewired sensor 22 mm x 22 mm (rectangular)	XN.TSP	20 pouches
	Infants and Children	Sensor with clip 30 mm (round)	XN.TSN	25 pouches
	Adults	Sensor with clip 45 mm (round)	XN.TSA	25 pouches

Products and options may not be available in every region or country. Please contact [mail@osypkamed.com](mailto:mail@osypkamed.com) for availability in your region and country and for further information.

### Osypka Medical GmbH

Albert-Einstein-Strasse 3  
12489 Berlin, Germany  
Phone: +49 (30) 6392 8300  
Fax: +49 (30) 6392 8301  
E-Mail: [mail@osypkamed.com](mailto:mail@osypkamed.com)  
[www.osypkamed.com](http://www.osypkamed.com)

### United States of America: Osypka Medical, Inc.

dba. Cardiotronic  
La Jolla, CA 92037, USA  
Phone: +1 (858) 454 0600  
Fax: +1 (858) 454 0640  
E-Mail: [mail@cardiotronic.net](mailto:mail@cardiotronic.net)  
[www.cardiotronic.net](http://www.cardiotronic.net)

U.S. Patent Nr. 6,511,438. Other patents pending.

iSense, AESCULON, ICON, ICONCore, Cardiotronic, Electrical Cardiometry, Window to the Circulation and logos are trademarks of Osypka Medical.

Copyright © 2024 Osypka Medical

