

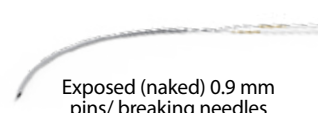
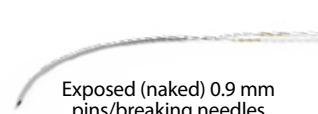
StimConnect™

for PACE® Model 300 with Unipolar Terminals

Extension Cables for Surgical Pacing Applications

PACE® Model 300 is primarily used for temporary surgical cardiac pacing applications. Surgical wires (temporary myocardial electrodes or heart wires) for temporary pacing are either equipped with 2 mm exposed (naked) plug pins or 0.9 mm pins/ needles. Alligator clips facilitate extraordinary temporary pacing applications.

Disposable cables are for single use, reusable extension cables can be sterilized 25 times by ethylene oxide or steam. A punch card allows the identification of the sterilization cycles.

Pacemaker	Disposable Extension Cable StimConnect™		Surgical Pacing Electrode
	Connection Pacemaker- Side	Connection Patient-Side	
	 Article No: F7818/BOSYLC (blue) F7818/OSYLC (white) Plugs with touch protected 2 mm pins fitting unipolar terminals of PACE® Model 300 Connects exposed (naked) pins with diameter between 0.9 to 2.0 mm Packaging Unit: 10 Pieces		 Adapters with exposed (naked) 2 mm pins  Exposed (naked) 0.9 mm pins/ breaking needles
	 Article No: F7833 Plugs with exposed 2 mm pins fitting unipolar terminals of PACE® Model 300 Covered alligator clips allow versatile connection Packaging Unit: 10 Pieces		
	Reusable Extension Cable StimConnect™  Article No: F7816/B33 (blue) F7816/33 (white) Plugs with exposed 2 mm pins fitting unipolar terminals of PACE® Model 300 Connects exposed (naked) pins with diameter between 0 to 2.0 mm Packaging Unit: 5 Pieces		 Adapters with exposed (naked) 2 mm pins  Exposed (naked) 0.9 mm pins/breaking needles

PACE® Model 300 with unipolar terminals accept- ing pins of 0.9 to 2.0 mm

(order number: P/N 2005)



Products and options may not be available or may vary by region and country. Please contact mail@osypkamed.com for availability in your region and country, for more information and for other StimConnect® extension cable variants.